



TabTimer.
helping to keep people on time

www.TabTimer.com.au

1300 TAB TIMER (1300 822 846)

APR 99 127 015 940

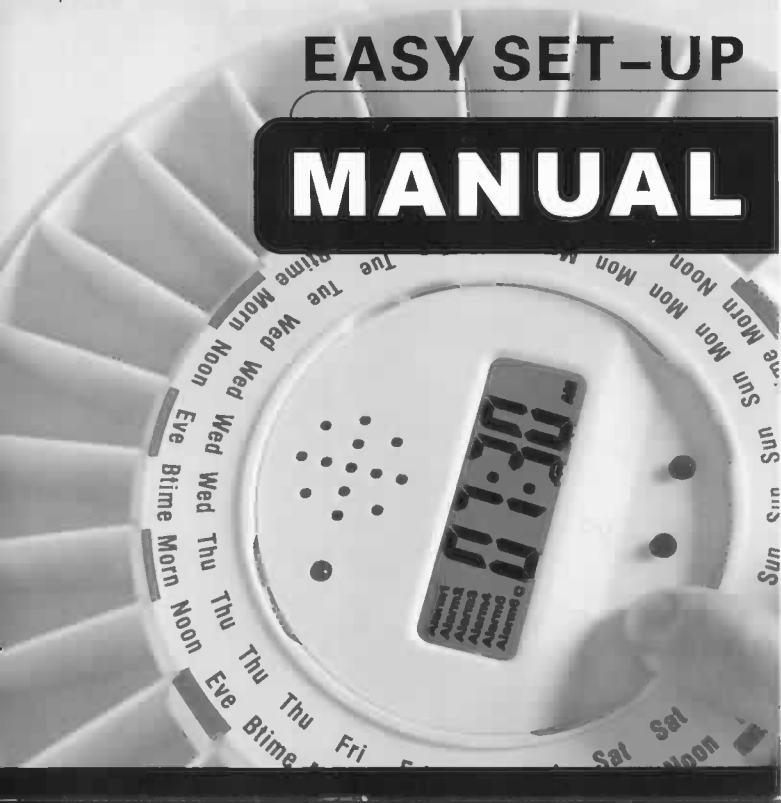


**AUTOMATIC
PILL DISPENSER**

Do not place Automatic Pill Dispenser
in direct sunlight or strong ultraviolet
light as it may cause malfunction
in operation.


EASY SET-UP

MANUAL





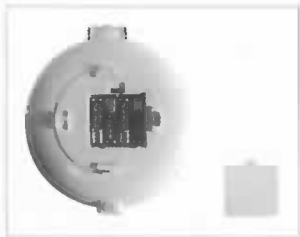1. Insert the batteries 2. How to open



Insert the batteries (4 - AA) into the battery compartment correctly.



Press gently with thumb on front tab of Med-E-Lert as shown.

3. How to clean the tray

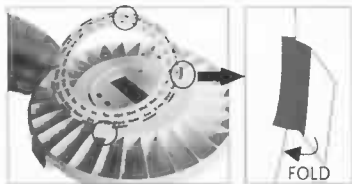
Clean only with damp cloth or sponge.



4. Loading the Med-E-Lert Dispenser



1. One time per day
2. Two times per day
3. Three times per day
4. Four times per day
5. Five times per day
6. Six times per day





The Med-E-Lert is supplied with 6 discs. The discs identify the day of the week, and the appropriate number of doses per day (1, 2, 3, 4, 5 or 6 doses per day).

The appropriate disc for the medication regime is then placed on the tray; fold tabs down to secure in pill compartment, prior to filling it, and used as a guide as to where to place the pills.

Organizing the medication, assigning dose times and determining what medications go into each slot is the responsibility of you, your caregiver, and/or your health care provider.

The refill frequency will depend upon the number of daily medications for example:

One time a day medication refill every 28 days

Twice a day medication refill every 14 days

Three times a day medication refill every 9 days

Four times a day medication refill every 7 days

Five times a day medication refill every 5 days

Six times a day medication refill every 4 days

5. 12/24-Hour Time Formats



MED-E-LERT automatically defaults to a 24 hour time format. To switch to a 12 hour format, press buttons 2 and 3 simultaneously. The AM or PM symbol will appear on the right side of the display when the clock is in the 12 hour time format. Release the buttons when the desired time format is displayed.

▲ Press button **2** and button **3** simultaneously



6. Setting the clock

(button 1- left; button 2- middle; button 3- right)

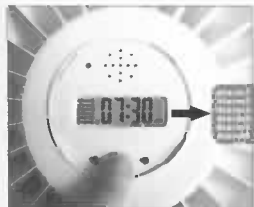


▲Press button ①

1. Press and hold button 1 for a few seconds (until the symbol ":" blinks) and then release.



BLINK




▲Press button ②

2. Set the correct time by pressing button 2 (hours) and button 3 (minutes).



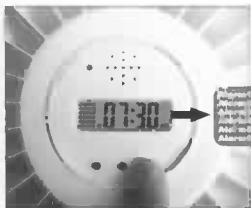BLINK
Set hour



▲Press button ③

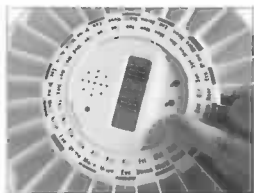
3. Store clock time by pressing button 1.




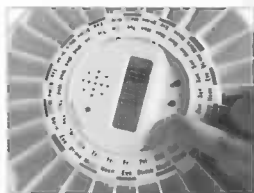
BLINK
Set minute




7. Setting and Checking the Alarm Times




1. Press button 1 and then release. The alarm 1 symbol will begin to blink. To set alarm time 1, press buttons 2 (hours) and 3 (minutes). If you wish to cancel the alarm, press button 2 until is displayed, which occurs after 23:00 or 11:00 p.m. 



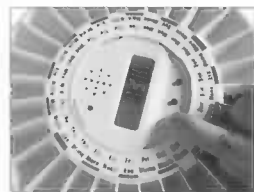2. After you have programmed the correct alarm time, press button 1 again to store the alarm time and to advance to the next Alarm setting.



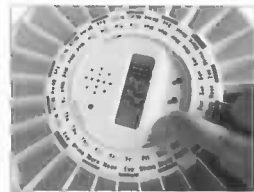3. The Alarm 2 symbol will blink. To set alarm time 2, use buttons 2 & 3.



4. Continue as above to set Alarm setting 3, 4, 5 and 6.

5. End by pressing button 1 until the Alarm symbols cease to blink.

- If buttons 2 or 3 are depressed, time display progresses automatically. If you let the display continue to blink in the setting mode, it will automatically return to normal function after one minute.

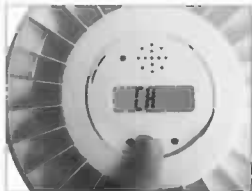
8. Alarm Tone(s)



▲ Press button 3 first, then press button 1 to choose the next alarm tone.

MED-E-LERT has 3 different alarm tones and also silence to accommodate different hearing levels. Symbols A, B, and C on the display indicate the selected alarm tone, if it be selected with "silence", then neither A nor B or C appears on the display. To change the alarm tone, first hold down button 3 and then, at same time press and hold down button 1 to choose the next alarm tone. Continue to hold button 3 and press button 1 to listen to and select the different tones.

9. Testing



▲ Press and hold button 2 for about 5 seconds, it display CH

MED-E-LERT has a built-in function to test the device. To enter the test mode, press button 2 for about 5 seconds, the CH should appear on the display. To exit test mode, press button 3.

CH Test Mode -When the display shows CH, pill box is ready to be tested. To perform a system function test, press button 1. The unit will proceed through the following system test modes.

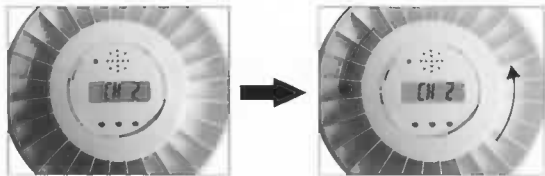
■ CH 1 Display Test

The display shows CH 1. All segments on the display shall light up.




■ CH 2 Display Test

The display shows CH 2. The motor shall rotate the tray one compartment.



■ CH 3 Alarm Signal Test

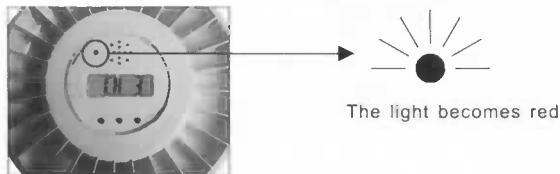The display shows CH 3. The alarm should sound and LED light lights up. Turn the pill box over to test RESET function (as to dispense medicine). The sound and LED light should be stopped. The CH 4 message will appear on the display when the pill box is not turned over within 10 seconds.



- If **MED-E-LERT** does not rotate properly the motor will stop after a few seconds and display the error message "CH6". Exit the test mode by pressing button 3.
- If there is a malfunction, please replace the batteries to reset.

10. How to Lock and unlock Med-E-Lert

Insert key and turn to left 1/4 turn to lock, Insert key and turn to right 1/4 turn to unlock.

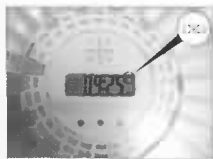


11. How to Stop Alarm



To stop the alarm and dispense the pills, pick up Med-E-Lert and tilt over, dumping the pills into the hand or a suitable container.

12. Low Battery Alert & Standby Wakeup



When the Low Battery indicator (as picture) appears on the LCD, the batteries are low and need replacing.

To save power, Med-E-Lert will stand by and shut display if no operation in several minutes. To wake it up, tilt it up or press the button 3.



TabTimer.
helping to keep people on time

200

www.TabTimer.com.au

1300 TAB TIMER (1300 822 846)

upply.com